



LAWRENCE
LIVERMORE
NATIONAL
LABORATORY

Probing Structure Property Relationships in Complex Engineering Silicones by ^1H NMR

S. C. Chinn, E. L. Gjersing, R. S. Maxwell, E.
Eastwood, D. Bowen, T. Stephens

October 2, 2006

American Chemical Society National Meeting
San Francisco, CA, United States
September 10, 2006 through September 14, 2006

Disclaimer

This document was prepared as an account of work sponsored by an agency of the United States Government. Neither the United States Government nor the University of California nor any of their employees, makes any warranty, express or implied, or assumes any legal liability or responsibility for the accuracy, completeness, or usefulness of any information, apparatus, product, or process disclosed, or represents that its use would not infringe privately owned rights. Reference herein to any specific commercial product, process, or service by trade name, trademark, manufacturer, or otherwise, does not necessarily constitute or imply its endorsement, recommendation, or favoring by the United States Government or the University of California. The views and opinions of authors expressed herein do not necessarily state or reflect those of the United States Government or the University of California, and shall not be used for advertising or product endorsement purposes.

PROBING STRUCTURE PROPERTY RELATIONSHIPS IN COMPLEX ENGINEERING SILICONES BY ^1H NMR

Sarah C. Chinn,¹ Erica Gjersing,¹ Robert S. Maxwell,¹ Eric Eastwood,² Dan Bowen,² Tom Stephens³

¹Lawrence Livermore National Laboratory, Livermore, CA 94550

²Honeywell FM&T Kansas City Plant, Kansas City, MO 64131

³Los Alamos National Laboratory, Los Alamos, NM, 87545

Introduction

It is generally accepted that the properties of polymeric materials are controlled by the network structure¹ and the reactions by which they have been constructed. These properties include the bulk moduli at creation, but also the properties as a function of age during use. In order to interpret mechanical properties and predict the time dependent changes in these properties, detailed knowledge of the effect of structural changes must be obtained. The degree and type of crosslinking, the molecular weight between crosslinks, the number of elastically ineffective chains (loops, dangling chain ends, sol-fraction) must be characterized.

A number of theoretical and experimental efforts have been reported in the last few years on model networks prepared by endlinking reactions² and the relationships of those structures with the ultimate mechanical properties.³ A range of experimental methods have been used to investigate structure including rheometric, scattering, infrared,²⁹Si MAS and CPMAS, ^1H relaxation measurements, and recently ^1H multiple quantum methods.

Characterization of the growth of multiple quantum coherences have recently been shown to provide detailed insight into silicone network structure by the ability to selective probe the individual components of the polymer network, such as the polymer-filler interface or network chains. We have employed recently developed MQ methods to investigate the structure-property relationships in a series of complex, endlinked filled-PDMS blends. Here, a systematic study of the relationship between the molecular formulation, as dictated by the amount and type of crosslinks present and by the remaining network chains, and the segmental dynamics as observed by MQ NMR was performed.

Experimental

Materials. The silicone-based materials were constructed using a series of varying molecular weight polydimethylsiloxane (PDMS) polymers, which were synthesized using a tin(II) ethylhexanoate catalyst and diatomaceous earth filler. Three different crosslinking agents were used: (1) polymethylhydrosiloxane (PMHS), which contains many (up to 60) crosslinking sites; (2) tetra-n-propoxysilane (TPS), containing up to four crosslinking sites; and (3) diphenylmethyl silanol (DPMS), which forms terminating crosslinks. The relative amounts of the various molecular weight PDMS components and the PMHS crosslinking agents were systematically varied, resulting in a series of formulations with different mechanical properties. An example of a final configuration is shown in Figure 1, which demonstrates the different types of crosslinks formed using PMHS vs. TPS.

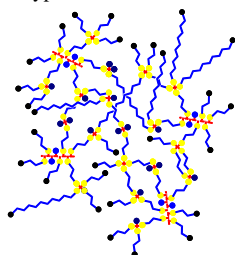


Figure 1. Typical structure of silicone elastomers showing different types of crosslinks and chain lengths.

^1H NMR Methods. Experiments were performed at 400.13 MHz on a Bruker Avance 400 spectrometer using a Bruker TBI (HCX) 5mm probe. In all cases, small (0.1 cm x 0.1 cm x 0.1cm) squares of elastomer were cut from a larger piece and set in the portion of a 5 mm NMR tube that would be within the coil volume of the probe. Traditional ^1H spin echo experiments were performed using 90° pulse lengths of 5.5 μs and recycle delay times of 6 sec, and echo delays were varied from 0-15 ms. Echo intensity curves were plotted on a logarithmic scale and fit to a bi-exponential decay representing two separate relaxation processes.

Multiple quantum NMR experiments were performed using the refocused multiple quantum excitation and reversion pulse sequence described previously.^{4,5} Pulse lengths of 5.5 μs were used with delay $\Delta 1$ and $\Delta 2$ equal to 4.83 μs and 6.16 μs respectively. Detailed descriptions of the experimental procedure can be found elsewhere.⁵

Results and Discussion

^1H relaxation and multiple quantum NMR experiments were first performed on a series of polymers prepared with varying PMHS concentrations. The spin-echo relaxation data was fit to an exponential decay to determine the T_2 relaxation time. As seen in Figure 2, the value for $1/T_2$, which can be correlated to the relative crosslink density or stiffness of the material,⁶ increased with increasing PMHS. MQ-NMR buildup curves were fit using two different values of the residual dipolar coupling (RDC).⁵ The higher RDC values represent regions of higher order, such as chains at a polymer-filler interface, while the lower RDC values result from chains with higher mobility, such as longer network chains. The ratio of RDC(low):RDC(high) can be used to provide information about the molecular partitioning of network vs. interfacial regions in a polymer system. In the case of varying PMHS, the ratio of RDC (low):RDC (high) decreased with increasing PMHS, indicating that higher concentrations of PMHS lead to larger contributions from the stiffer portion of the polymer. The trends observed in the NMR data follow closely with mechanical testing of load requirements, indicating that the amount of load increases with increasing PMHS concentration.

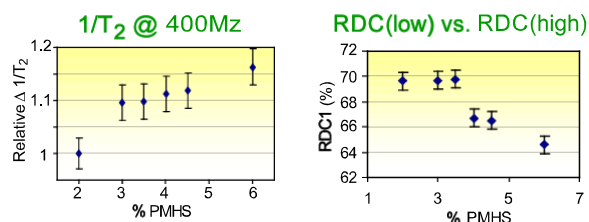


Figure 2. NMR results of variable PMHS study.

Comparable studies were performed with samples of varying molecular weights of PMHS crosslinkers. In general, a slight increase in $1/T_2$ was observed at higher molecular weights, as the number of available crosslinking sites increased. Additionally, both the low and high contributions to the RDC's were observed to decrease with higher MW of the PMHS, as did the ratio of RDC(low):RDC(high). Again, these changes are due to the presence of more crosslinking sites with higher MW PMHS chains.

Detailed investigations were also performed on hybrid materials consisting of varying amounts of low-MW, high-MW, or mid-MW (termed "low-high" MW) PDMS components. In these studies, the relative amounts of high-MW and low-high-MW PDMS were varied. NMR relaxation experiments revealed that the value of $1/T_2$ decreased with increasing ratios of H/LH MW materials due to a softening of the material. While there were few significant changes in the high or low contributions to the residual dipolar coupling, the ratio of RDC(low):RDC(high) increased with increasing ratios of H/LW MW PDMS. It is speculated that the relative populations of the mobile regions, which contribute to RDC(low) and also ultimately to the T_2 value, generally increase with increasing numbers of long PDMS chains. The presence of the longer chains leads to a softer polymer. Again, this is verified by mechanical load measurements, which indicate that the load decreases with increasing ratio of H/LH MW PDMS.

Conclusions

The material properties of silicone materials can be tailored by varying proportions of individual components of the polymer. ^1H relaxation and multiple-quantum NMR can be used to reveal network-specific information and provides valuable information about the segmental dynamics as a function of formulation. In the study presented here, increasing amounts of the PMHS crosslinker leads to a stiffer polymer, while larger amounts of the high molecular weight PDMS chains contributes to a softening of the polymer. A detailed focus on chemistry followed by detailed characterization of the material properties of polymer formulations is crucial in achieving a thorough understanding of the structure-property relationships of the materials.

Acknowledgements. This work was performed under the auspices of the U.S. Department of Energy by UC, Lawrence Livermore National Laboratory under contract # W-7405-ENG-48.

References

- Larsen, A.L., Hansen, K., Sommer-Larsen, P., Hassager, O., Bach, A., Ndoni, S., Jørgensen, *Macromolecules* **2003**, *36*, 10063.
- Gilra, N., Cohen, C., Panagiotopoulos, *J. Chem. Phys.* **2000**, *112* 6910.
- Patel, S.K., Malone, S., Cohen, C., Gillmor, J.R., Colby, R., *Macromolecules* **1992**, *25*, 5241.
- Saalwächter, K., *J. Chem. Phys.* **2004**, *120*, 454.
- Chinn, S., DeTeresa, S., Sawvel, A., Shields, A., Balazs, B., Maxwell, R.S., *Poly. Degrad. Stab.*, **2006**, *91*, 555.
- Maxwell, R.S., Balazs, B., *J. Chem. Phys.* **2002**, *116*, 10492.